


# Steroids

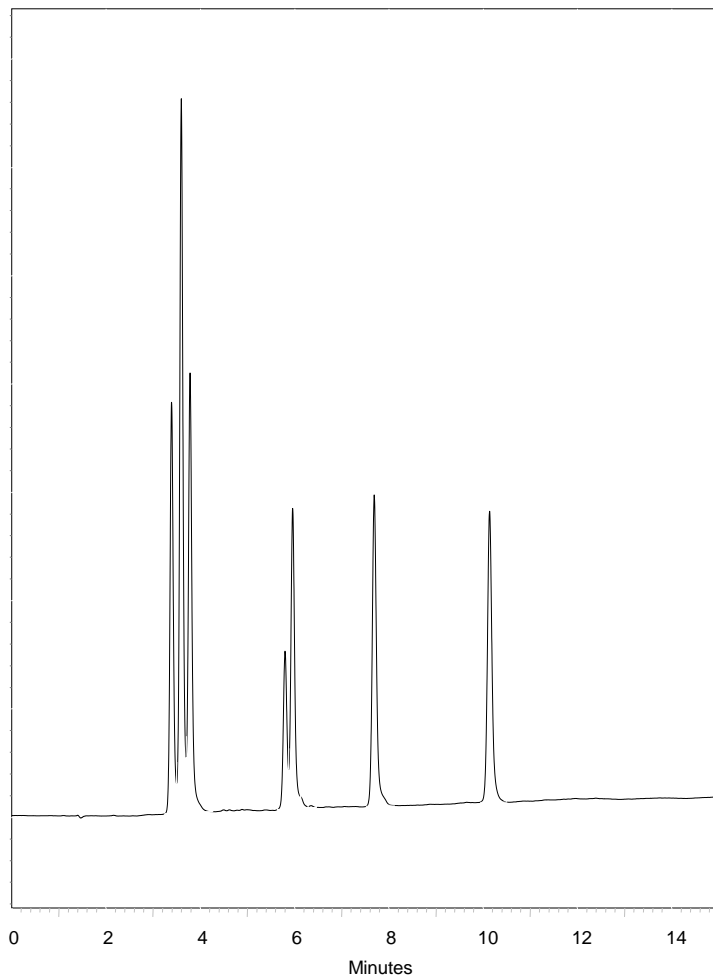

Column : **5 $\mu$  Fortis C18 100x4.6mm**  
p/n : **F18-050505**

Mobile Phase: **A: H<sub>2</sub>O + 0.1% Formic acid**  
**B: ACN + 0.1% Formic acid**

Gradient: **30-100% in 10 minutes**

Flow : **1.0ml/min**

Temp : **20°C**

$\lambda$  : **254nm**

- 1. Prednisolone**
- 2. Prednisone**
- 3. Cortisone**
- 4. 11 $\alpha$  Hydroxyprogesterone**
- 5. Hydroxytisone-21-acetate**
- 6. 17 $\alpha$  Hydroxyprogesterone**
- 7. Progesterone**