


Herbicides

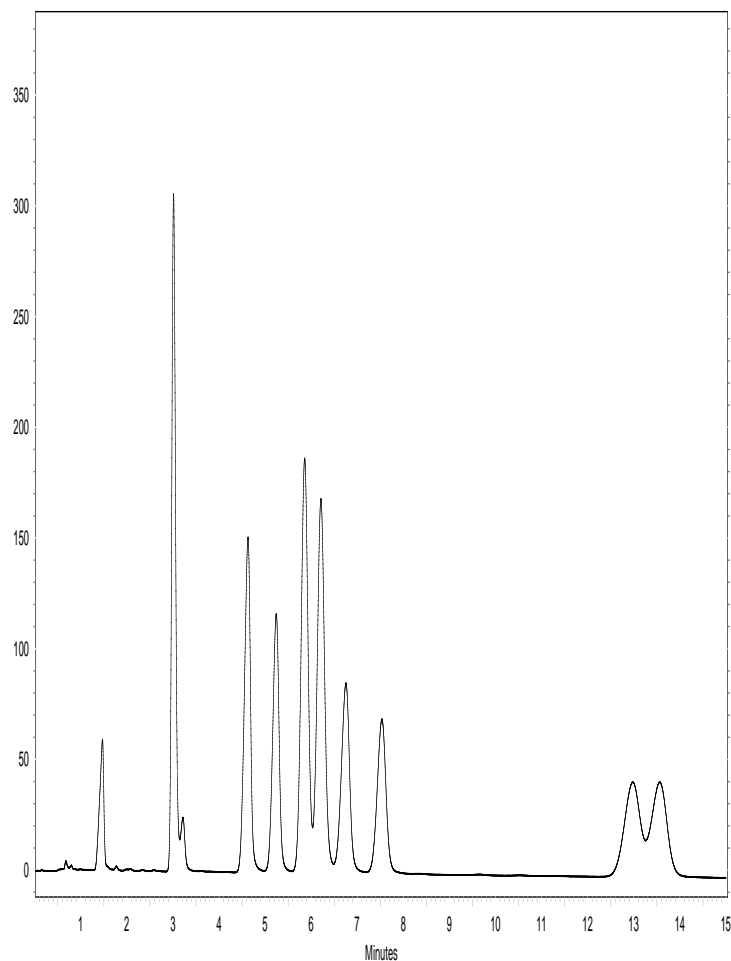

Column : **3 μ Fortis Cyano 50x2.1mm**
p/n : **FCN-020303**

Mobile Phase: **80:20 –**
H₂O + 0.2% Acetic acid:ACN

Flow : **0.2ml/min**
Temp : **25°C**
 λ : **280nm**

- 1. Banzel**
- 2. 2,4-D**
- 3. MCPA**
- 4. PCOC**
- 5. 2,4-DCP**
- 6. 2,4-DP**
- 7. CMPP**
- 8. 2,4-DB**
- 9. MCPB**