


# Sulfa Drugs

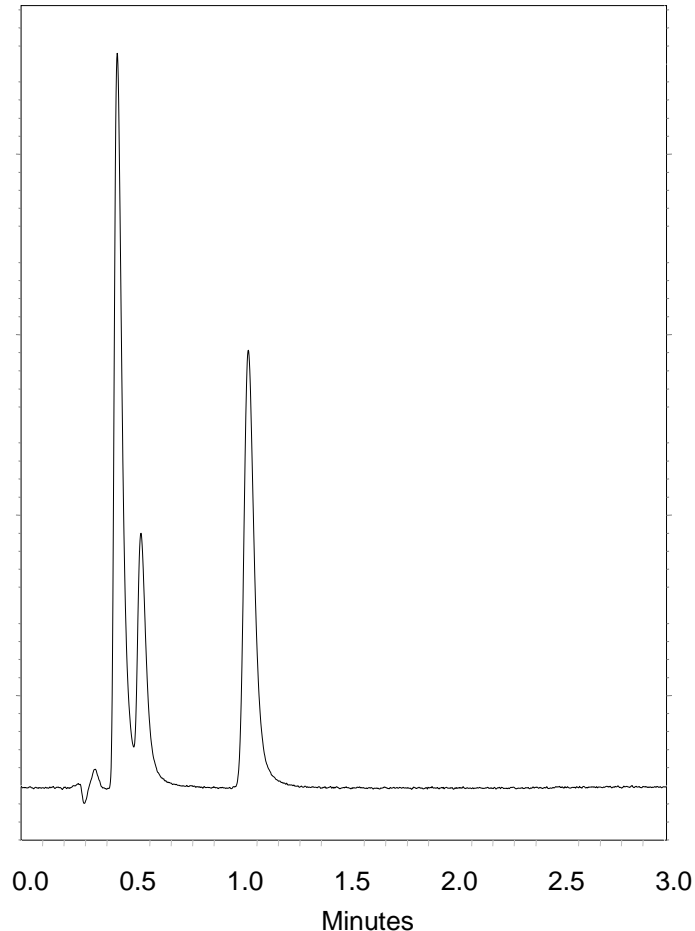

Column : **3 $\mu$  Fortis Pace H<sub>2</sub>O  
30x2.1mm**

p/n : **FHO-020203**

Mobile Phase: **70:30 H<sub>2</sub>O + 0.1% Formic acid  
:ACN + 0.1% Formic acid**

Flow : **0.2ml/min**

Temp : **25°C**

$\lambda$  : **254nm**

- 1. Sulfathiazole**
- 2. Sulfamerazine**
- 3. Sulfamethoxazole**