


Tyrosine Metabolism

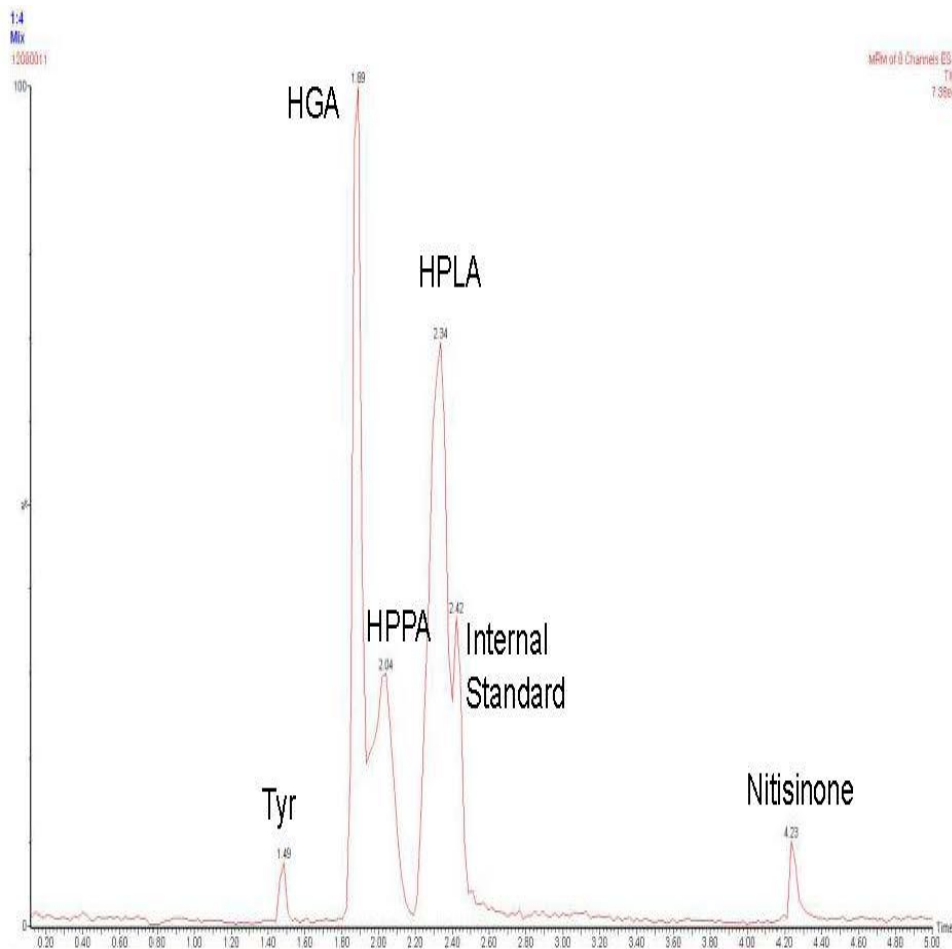

Column : **3 μ Fortis H₂O**

100x4.6mm

p/n : **FHO-050503**

Mobile Phase: **85:15 ACN:Water**

Flow : **1.0ml/min**

Temp : **20°C**

λ : **254nm**

1. Tyrosine
2. Homogentisic acid (HGA)
3. Benzoquinone acetic acid (BQA)
4. Hydroxyphenylpyruvic acid (HPPA)
5. Hydroxyphenyllactic acid (HPLA)
6. Nitrosinone