


# Substituted Benzenes



- Column : **Evosphere C18-PFP 150x4.6 5um**
  - p/n : **EVO18FP-050705**
  - Mobile Phase: **50:50 Water:MeOH**
  - Flow : **1.0 ml/min**
  - Temp : **20°C**
  - $\lambda$  : **210nm**
1. **1,2,3 Trimethoxybenzene**
  2. **1,2 Dimethoxybenzene**
  3. **1,2,4 Trimethoxybenzene**
  4. **1,4 Dimethoxybenzene**
  5. **Anisole**
  6. **1,3 Dimethoxybenzene**
  7. **1,3,5 Trimethoxybenzene**
  8. **Toluene**

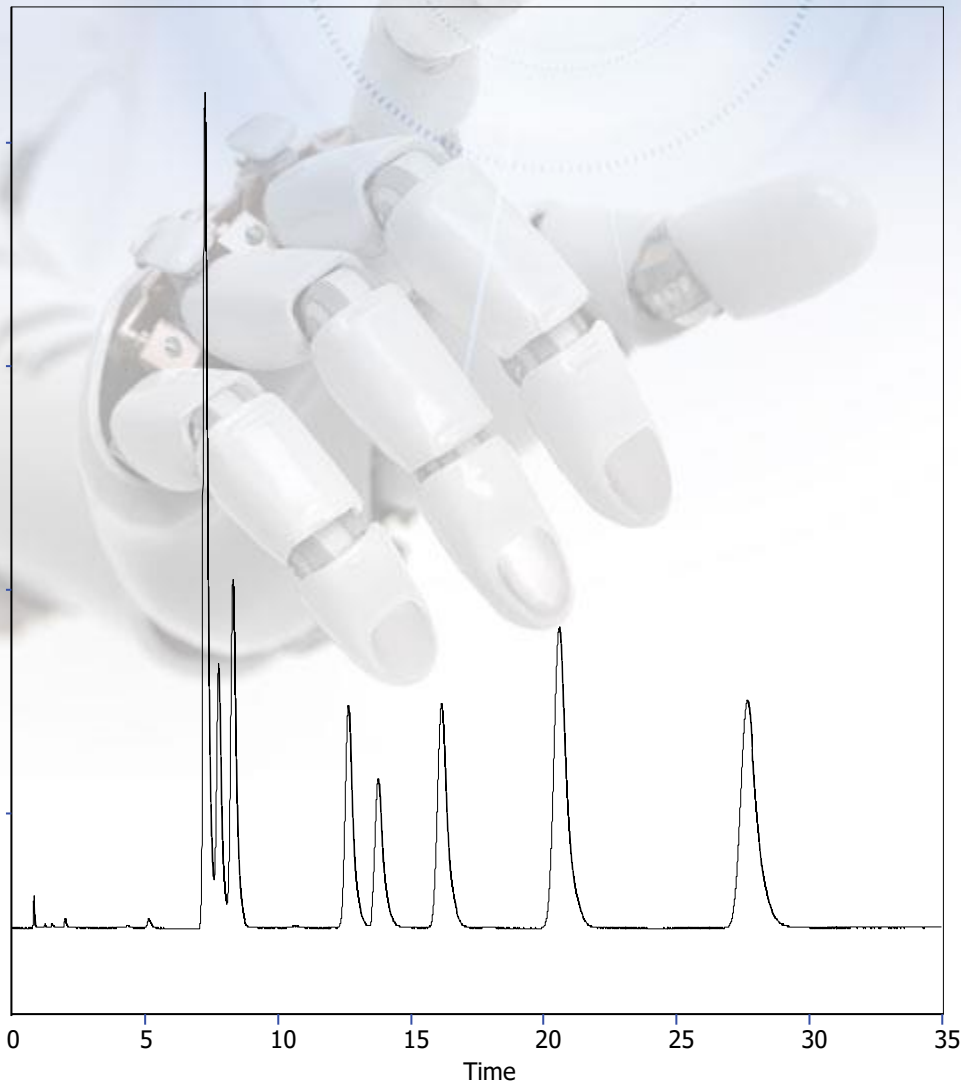