


Steroids



- Column : **Evosphere AQUA 150x4.6 3um**
 - p/n : **EVOAQ-050703**
 - Mobile Phase:
 - A- Water + 0.1% Formic acid**
 - B- ACN + 0.1% Formic acid**
 - 30-100% B in 10 mins**
 - Flow : **1.0 ml/min**
 - Temp : **AMB**
 - λ : **254nm**
1. **Prednisolone**
 2. **Prednisone**
 3. **Cortisone**
 4. **17 α Hydroxyprogesterone**

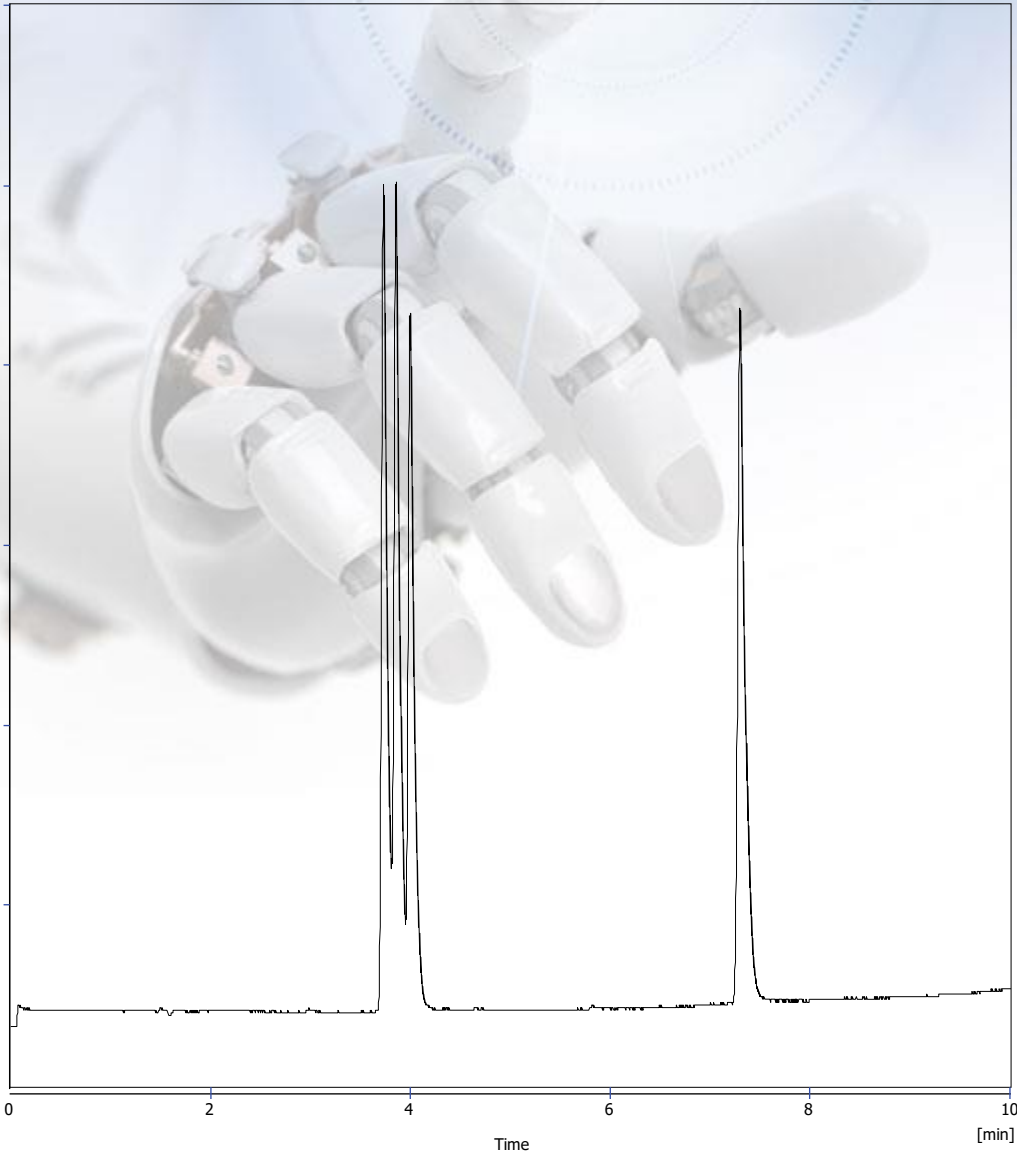