


# Polar Test Mix



- Column : **Evosphere AQUA 150x4.6 3um**
  - p/n : **EVOAQ-050703**
  - Mobile Phase:  
**75:25 20mM KH<sub>2</sub>PO<sub>4</sub>(pH7.6):ACN**
  - Flow : **1.0 ml/min**
  - Temp : **30°C**
  - $\lambda$  : **270nm**
1. **4-Nitrocatechol**
  2. **Resorcinol**
  3. **3-Methoxycatechol**
  4. **4-Methylcatechol**

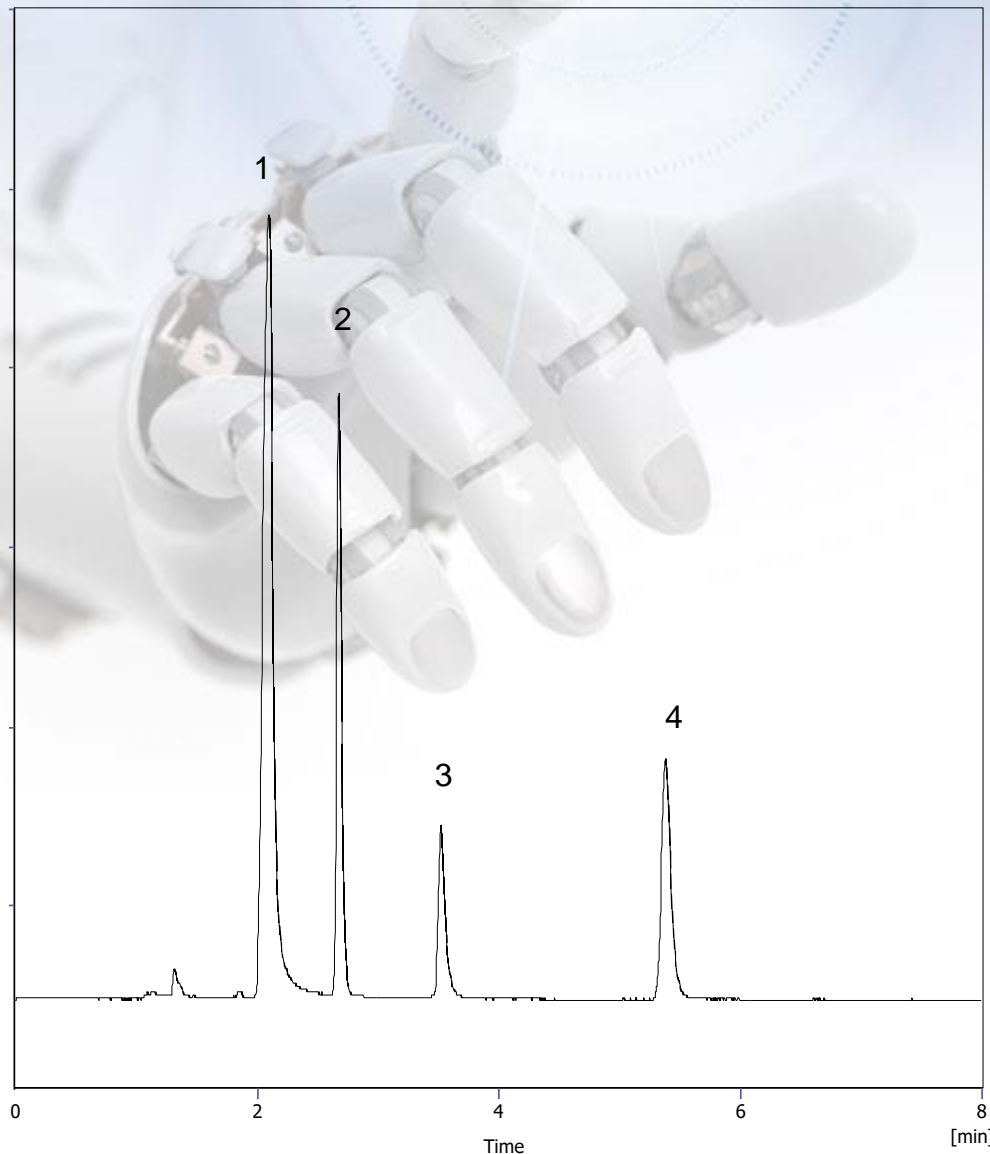