


# $\beta$ -Antagonists and Diuretics



- Column : **Evosphere C18/AR 100x2.1x3 $\mu$**
- p/n : **EVO18AR-020503**
- Mobile Phase:
- **A- 20mM Potassium Phosphate at pH 2.7**
- **B- 65:35 20mM Potassium Phosphate at pH 2.7: ACN**
- Gradient: **50%- 95% B in 5 minutes**
- Flow : **0.2 ml/min**
- Temp : **36°C**
- $\lambda$  : **230nm**

1. **Atenolol**
2. **Pindolol**
3. **Acetazolamide**
4. **Furosemide**
5. **Propranolol**
6. **Bumetanide**

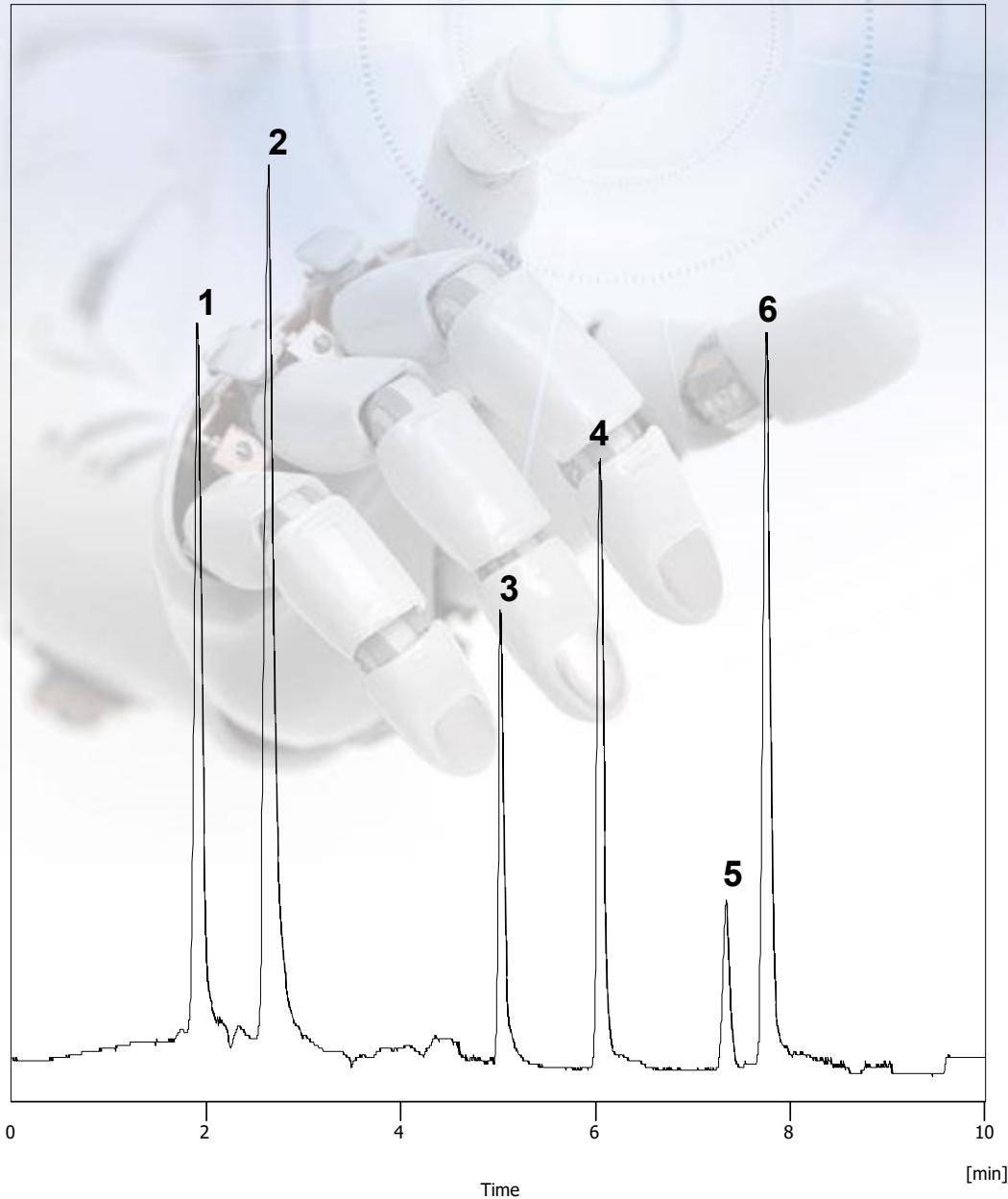