



KFDA 83 59

➤ Shimadzu Nexera with LCMS-2020

➤ Shimpack VP-ODS

250x4.6mm 5um

1) ?

2) ?

3) ?

4) ?

➤ A: 5mm ammonium acetate

➤ B: 5mm ammonium acetate in MeOH

Time	%B
------	----

0	20
---	----

0.3	20
-----	----

5.5	95
-----	----

7.5	95
-----	----

10	20
----	----

➤ 0.45ml/min


➤ 50oC

➤ 360nm




KFDA 83 59

Shimpack VP-ODS 150x2.1mm 2.2um

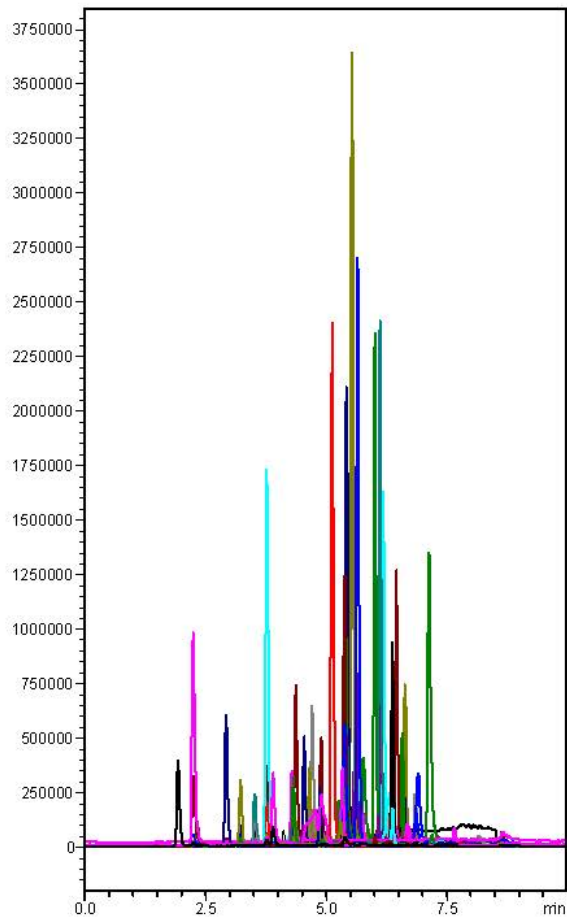

1.7um Fortis C18 50x2.1mm

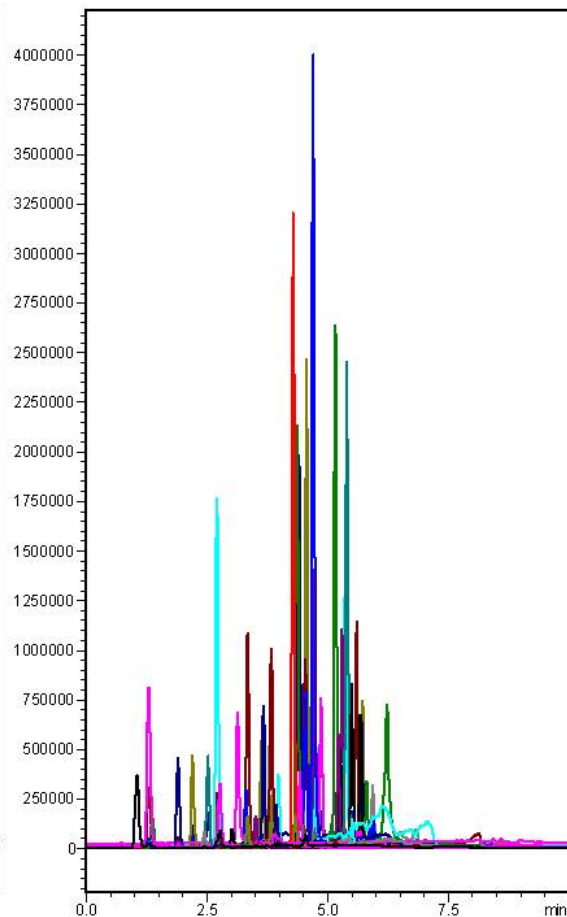

1.7um Fortis C18 30x2.1mm

