



PAH's

- Shimadzu LC-20A
- Shimpack VP-ODS 250x4.6mm 5um
- 1.7um Fortis C18 50x2.1mm

➤ A: Water

➤ B: ACN

Time	%B
0	70
20	70
27	100
30	70
40	70


- 1.0ml/min
- 40oC
- 254nm

- 1) Naphthalene
- 2) Acenaphthylene
- 3) Acenaphthalene
- 4) Fluorene
- 5) Phenanthrene
- 6) Anthracene
- 7) Fluoranthene
- 8) Pyrene
- 9) Benz[a]anthracene
- 10) Chrysene
- 11) Benzo[b]fluoranthene
- 12) Benzo[k]fluoranthene
- 13) Benzo[a]pyrene
- 14) Diben[a,b]anthracene
- 15) Benzo[ghi]perylene
- 16) Indeno[1,2,3,cd]pyrene



PAH's

Shimpack VP-ODS 250x4.6mm 5um

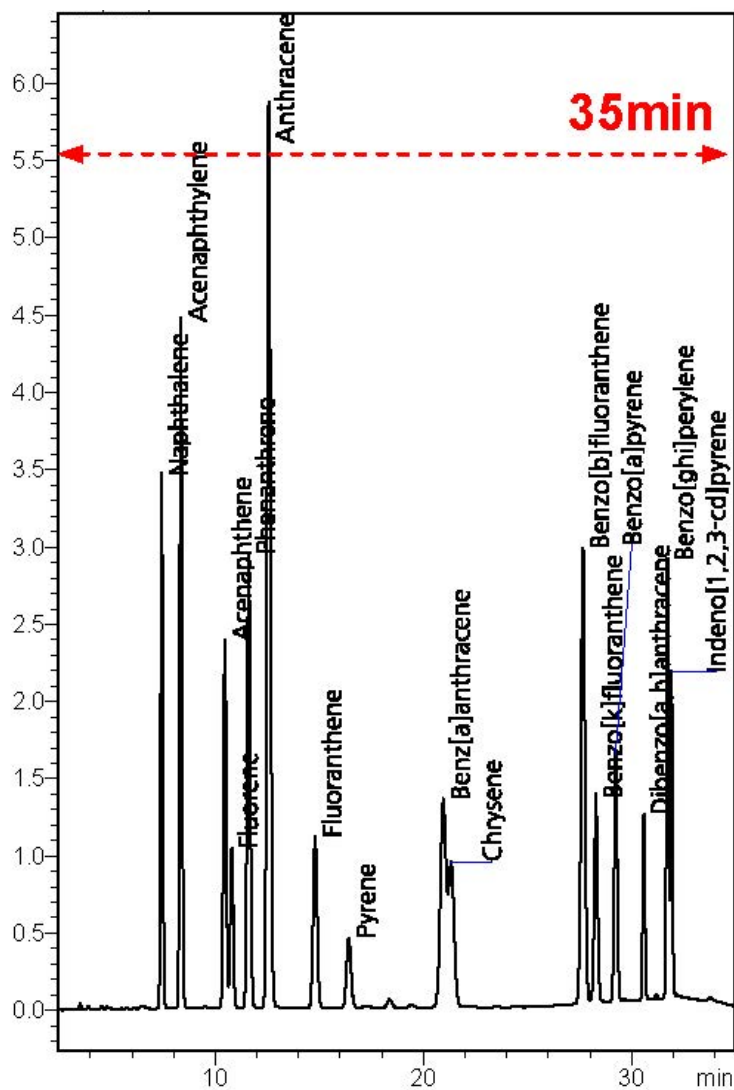

1.7um Fortis C18 50x2.1mm

