

SpeedCore Diphenyl


- Column : **SpeedCore Diphenyl 150x4.6mm 2.6 μ**
- p/n : **SPPH-050726**

- Mobile Phase:
40:60 Water : MeOH


- Flow : **1.2ml/min**
- Temp : **40°C**
- λ : **210nm**

1. **4-Hydroxyestradiol (mw=288.38)**

2. **2-Hydroxyestradiol (mw=288.38)**

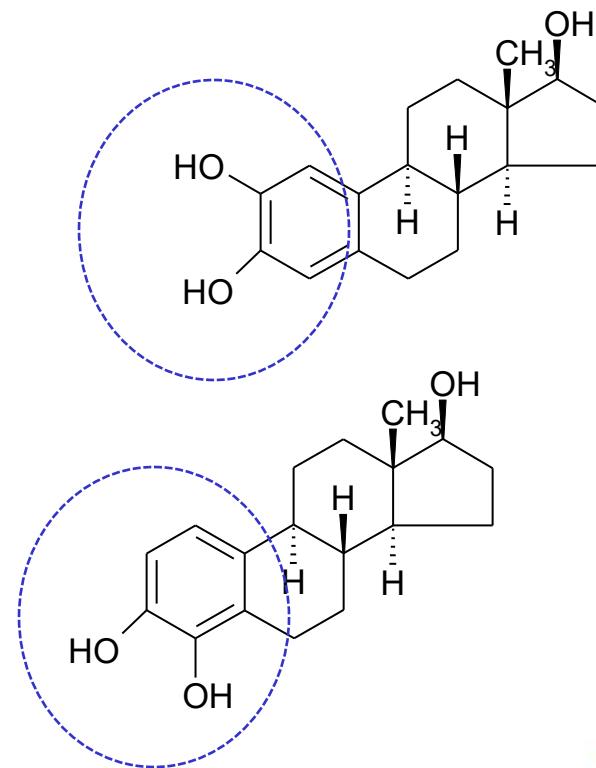

SpeedCore Diphenyl - Isomers

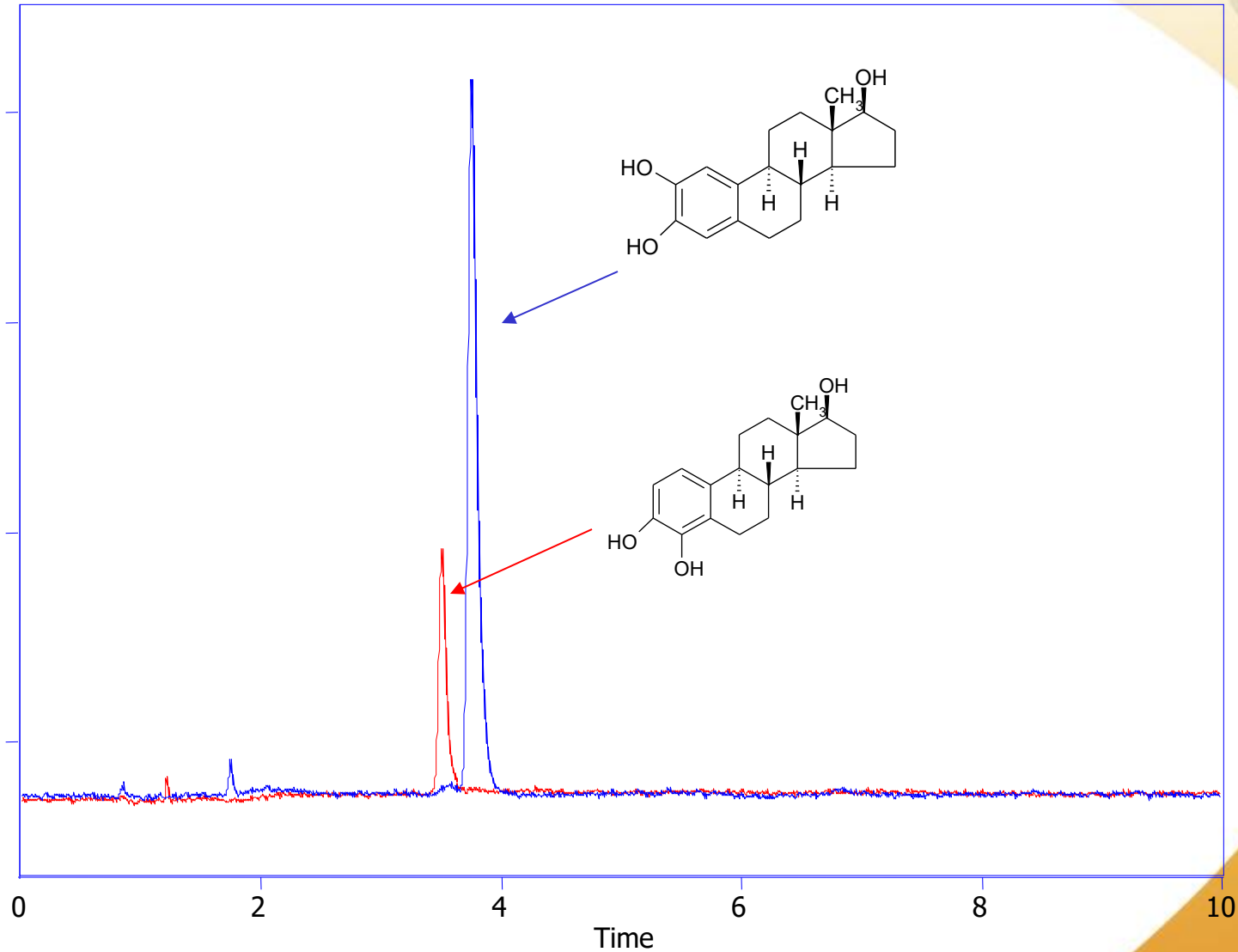

SpeedCore PFP - Isomers

