




Fortis BIO

Peptide Mix

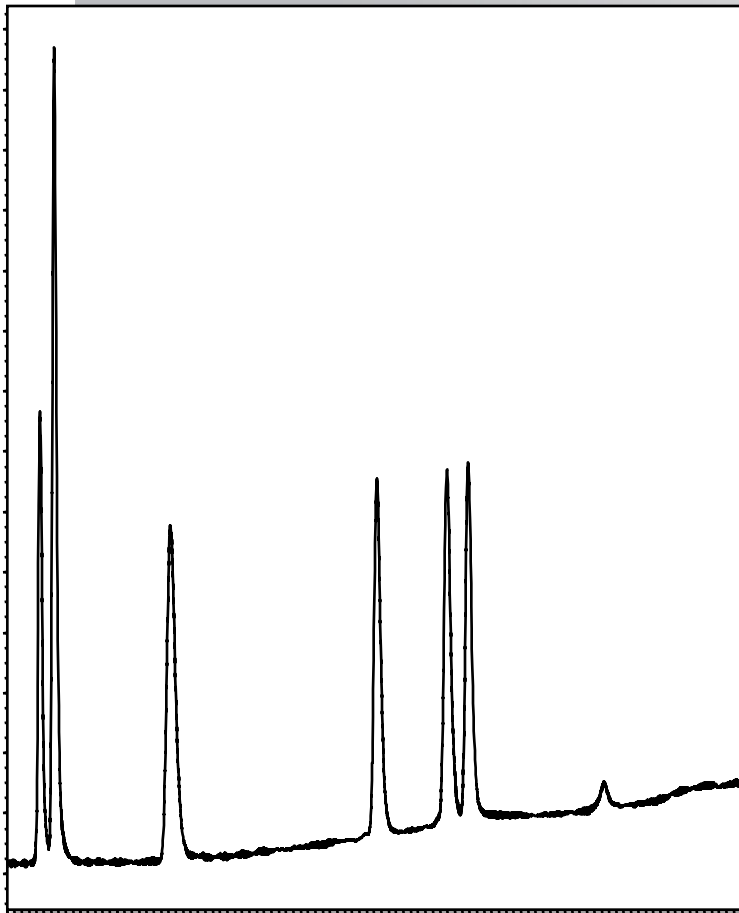

- Column : **FortisBioC18 50x2.1mm 5 μ**
- p/n : **BIO318-020205**
- Mobile Phase:
 - 10-40% B in 10mins**
 - A: 0.1% TFA**
 - B: ACN**
- Flow : **0.4 ml/min**
- Temp : **25°C**
- λ : **220nm**

1. **GLY-TYR (238.2 MW)**
2. **VAL-TYR-VAL (379.5 MW)**
3. **MET-Enkephalin (573.7 MW)**
4. **LEU –Enkephalin (555.6 MW)**
5. **Angiotensin II (1046.2 MW)**