




Fortis BIO

Protein Mix

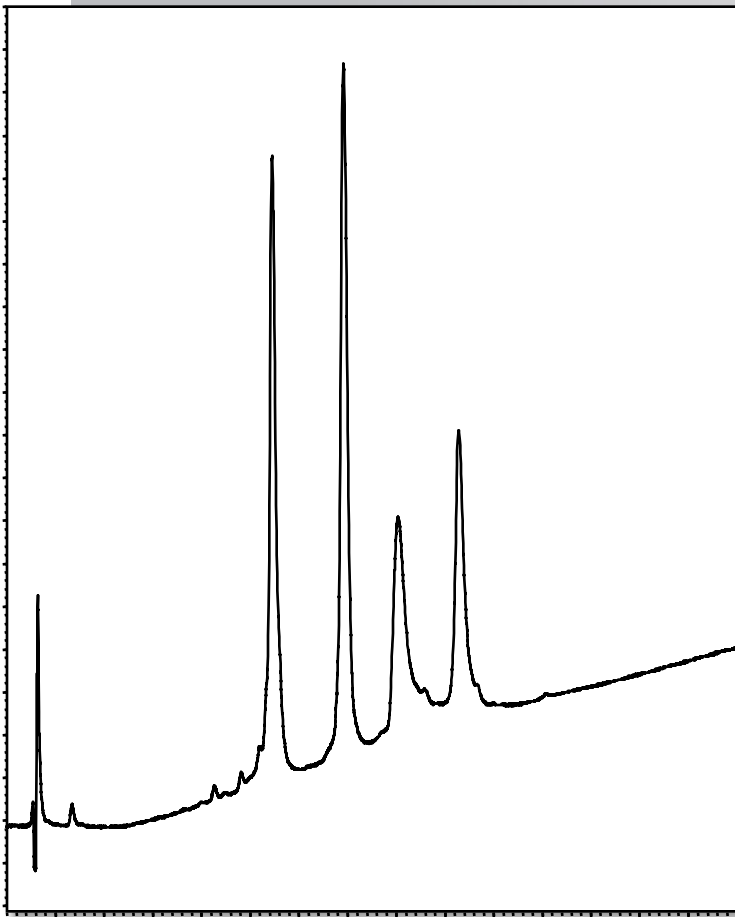

- Column : **FortisBioC18 50x2.1mm 5 μ**
- p/n : **BIO318-020205**
- Mobile Phase:
10-60% B in 15mins
A: 0.1% Formic acid
B: ACN
- Flow : **0.4 ml/min**
- Temp : **25°C**
- λ : **220nm**

- 1. Ribonuclease A**
- 2. Cytochrome C**
- 3. Holo- Transferrin**
- 4. Apomyoglobin**



Fortis™
Technologies Ltd.